




## MATERIAL PROCESSING

### LAB CORE CRUSHER

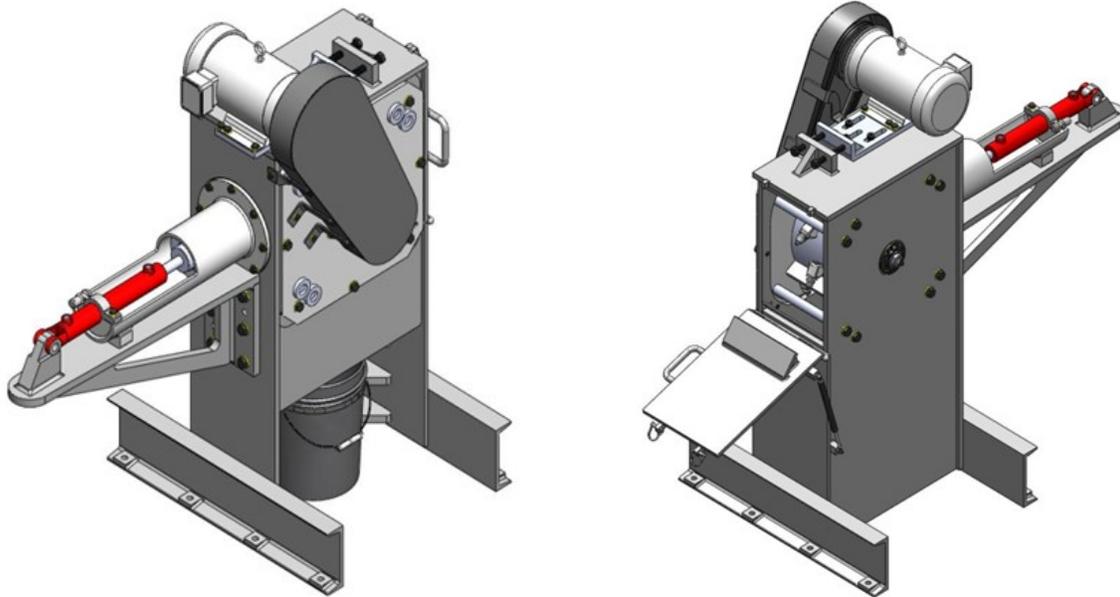

- Creates material that is representative of field milling
- The crusher's drum design replicates an actual milling machine's drum with similar wrap design, tooth spacing and drum RPM
- Saves time by crushing cores in seconds and eliminates the need for a Jaw crusher



Sieve size	As-crushed	Coarse gradation band
1"	99.1	85 – 100
¾"	94.4	75 – 92
No. 4	34.5	30 – 45
No. 30	6.3	1 – 7
No. 200 (not washed)	0.52	0.1 - 3